



**THE PRO PERSONAL THERAPY SYSTEM
PAIN RELIEF - MEDICATION FREE**



APPROVED MEDICAL PRODUCT
MANUFACTURED IN GERMANY EXCLUSIVELY FOR CYCLO-SSAGE

CE 0633

WWW.CYCLO-SSAGE.COM

THE PRO PERSONAL THERAPY SYSTEM

MEDICALLY
APPROVED



CYCLODIAL THERAPY
WITH
INFRARED HEAT

CHOICE OF THERAPY POSITIONS



Full Body Therapy Position



Foot Therapy Position



Prone Therapy Position

COMPONENT PARTS OF THE PRO PERSONAL THERAPY SYSTEM



DIGITAL HAND CONTROL



SHOULDER/UPPER BACK INFRARED
HEAT CUSHION



LOWER BACK INFRARED
HEAT CUSHION



FOOT/CALF INFRARED
HEAT CUSHION



12V CAR/CARAVAN
ADAPTOR




PORTABLE CARRY CASE



SMALL SUPPORT CUSHION



UPPER BODY ELEVATION CUSHION

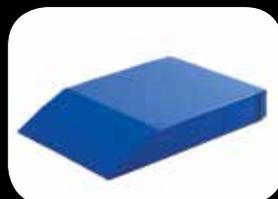

LEG/FOOT ELEVATION CUSHION




100/240V - 12V
MAINS TRANSFORMER

THERAPY PROGRAMME CHART & DIGITAL HAND CONTROL

Therapy programme chart showing the nine pre-selected therapy combinations that are available to select for your personal treatment.

CONDITIONS / SYMPTOMS	DISPLAY	15 MIN'S	30 MIN'S	45 MIN'S	60 MIN'S
THERAPY PROGRAMME 1					
tension and stiffness	1.1	2 min P1	4 min P1	6 min P1	8 min P1
general back ache	1.2	6 min P3	12 min P3	18 min P3	24 min P3
tired and aching muscles	1.3	7 min P2	14 min P2	21 min P2	28 min P2
THERAPY PROGRAMME 2					
sciatica	2.1	1 min P1	2 min P1	3 min P1	4 min P1
lower back pain	2.2	5 min P4	10 min P4	15 min P4	20 min P4
lumbago	2.3	5 min P2	10 min P2	15 min P2	20 min P2
muscle spasm	2.4	4 min P3	8 min P3	12 min P3	16 min P3
THERAPY PROGRAMME 3					
headaches and migraine	3.1	3 min P1	6 min P1	9 min P1	12 min P1
neck and shoulder tension	3.2	3 min P5	6 min P5	9 min P5	12 min P5
stress	3.3	9 min P3	18 min P3	27 min P3	36 min P3
THERAPY PROGRAMME 4					
arthritis and rheumatism	4.1	5 min P1	10 min P1	15 min P1	20 min P1
parkinson's and MS – fibromyalgia and ME	4.2	3 min P2	6 min P2	9 min P2	12 min P2
osteoporosis - Joint pain and muscle ache	4.3	7 min P3	14 min P3	21 min P3	28 min P3
THERAPY PROGRAMME 5					
sluggish lymphatic system	5.1	1 min P1	2 min P1	3 min P1	4 min P1
excess fluid	5.2	5 min P3	10 min P3	15 min P3	20 min P3
detoxification	5.3	4 min P2	8 min P2	12 min P2	16 min P2
cellulite	5.4	5 min P3	10 min P3	15 min P3	20 min P3
THERAPY PROGRAMME 6					
swollen legs and oedema	6.1	2 min P1	4 min P1	6 min P1	8 min P1
fluid retention	6.2	13 min P3	26 min P3	39 min P3	52 min P3
THERAPY PROGRAMME 7					
sport injuries	7.1	3 min P2	6 min P2	9 min P2	12 min P2
excess lactic acid and metabolic waste	7.2	3 min P4	6 min P4	9 min P4	12 min P4
strained and aching muscles	7.3	4 min P2	8 min P2	12 min P2	16 min P2
muscle spasm and cramp	7.4	5 min P3	10 min P3	15 min P3	20 min P3
THERAPY PROGRAMME 8 (to be used while lying on the stomach)					
digestive and abdominal problems	8.1	8 min P2	16 min P2	24 min P2	32 min P2
colitis and irritable bowel	8.2	2 min P4	4 min P4	6 min P4	8 min P4
period pains	8.3	5 min P2	10 min P2	15 min P2	20 min P2
THERAPY PROGRAMME 9					
breathing and respiratory problems	9.1	1 min P1	2 min P1	3 min P1	4 min P1
stress and anxiety	9.2	7 min P3	14 min P3	21 min P3	28 min P3
insomnia and sleep disorders	9.3	7 min P5	14 min P5	21 min P5	28 min P5

The beneficial effects of the therapy are cumulative – the more often these treatments are applied the better the results and the longer the effects will last.

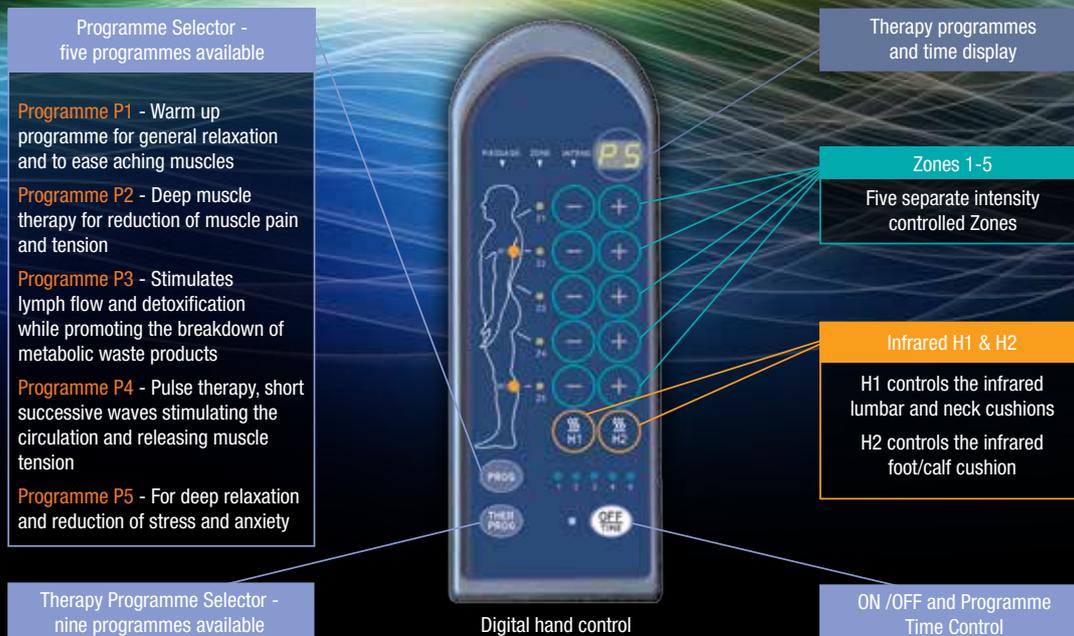

TECHNICAL DATA

Primary voltage: 100-240V AC 50/60 HZ
Secondary voltage: 12v DC / 8.33 Amps
Power consumption: Approx - 85W

Dimensions unfolded: L-195cm.W-60cm
Total Postal weight inc 2 boxes: Approx - 16.5 kg
Total Nett weight inc - accessories: Approx - 13.5 kg

SOME OF OUR MANY HAPPY CUSTOMERS



Mrs Tipping
Migraine and ME



Mr Hargreaves
Parkinson's Disease



Mr Headford
Rheumatoid Arthritis



Ms Speight
Scoliosis



Mrs Norris
Sciatica



Mrs Ansen
Osteoarthritis



Mrs Beumont
Plantar Fasciitis



Mrs Hutchinson
Fibromyalgia



Miss O'Connor
Lupus



Mrs Sisodia
Oedema



Mr Grant
Leg Amputation



Mr Luff
Broken Shoulder

SOME OF OUR MANY CUSTOMER TESTIMONIALS

Arthritis - *Anne Weaver*

"Just a quick email to say thank you unlike most products yours has done exactly what you told us it would to say that it has made a difference to my life is an understatement.....I have ceased taking pain killers as I don't need them..."

ME - *Kate Louise Bowskill*

"I have had ME/Fibro for 21 years and have been using my massage bed for a year and half now. During that time I have come off all pain medication. I was devastated when I found out I couldn't use it whilst pregnant but after seeking advice from my GP I am using it again and extremely glad I sought his advice."

FM (also Scoliosis, Spinal Fusion and ME) - *Mary Pooley*

"...The Cyclo-ssage has provided continued relief and comfort since I purchased it and is sometimes the only place I can go for relief. I love that there are no side effects (like with medications), that it is there at home when ever I need it and the quality of the product which still looks as good as new over two years on."

Back and Hip Pain - *Kenneth Alan*

"I purchased my Cyclo-ssage Pro Personal Therapy System from your company about two years ago and I would like to express how pleased I am with it.....it facilitates my movement and the relief from back and hip pain is wonderful. I know that it will not cure my medical condition but it certainly makes my life more comfortable. I cannot imagine life without it!"

Degenerative Facet Joint Disease - *Mike Waters*

"I have not used my walking stick since using the bed, my depression has lifted, my memory improved, I am sleeping better and I have given up smoking as well. In short the improvement that I have had since the massage bed is nothing short of miraculous."

TO VIEW MORE GO TO WWW.CYCLO-SSAGE.COM



Lady Jackson
Broken Back



Mr Barclay
Diabetes



Mrs Downs
Polymyalgia



Mrs Hart
High Blood Pressure



Mrs Mortimer
Spinal Curvature



Mr Reynolds
Broken Hip



Lewis
Cerebral Palsy



Mr Guy
Muscle Spasm



Mr Frodsham
Spinal Stenosis



Mr Elder
Spinal Condition



Mrs Brooks Stockton
Osteoarthritis - Neck



Mrs Tapley
Peripheral Neuropathy

TO VIEW MORE GO TO WWW.CYCLO-SSAGE.COM

Neck and Shoulder Pain - *Doreen Davenport*

"My husband suffers from arthritis and has a lot of pain in his neck and shoulders and the system helps relieve this pain. It is easy to use and I think it is absolutely marvellous!"

MS - *Susan Shand*

"...This morning I put a load of washing in the machine and managed to stand up again without my knees locking and me having to shout for some one to stand me up! I feel more supple and I also have a bit more feeling in my fingers...it's such a relaxing therapy and when the massage is finished you feel this wave going through your body, it's fab."

Lymphodema (also Poor Circulation and Back Pain) - *Caroline Waring*

"...Being in a wheel chair meant that my ankles used to swell because the circulation couldn't go back up my legs, but by using the System on the programme that helps lymphatic drainage that has now gone as well. It's very popular in the house and we almost need an appointment book for it now – I wouldn't hesitate to recommend it."

Low Energy Levels (also Neck and Shoulder Pain) - *Dianne Hughes*

"...I feel so much more free and comfortable, particularly in my lower back, neck and shoulders. I have also noticed a difference in how I feel at work as well. My husband was shocked at first when I told him I was buying it but now he uses it on a regular basis and I have even found him asleep on it before!"

Parkinson's Disease - *Stella Dowding*

"Results are good, my walking has improved a great deal, as have all my movements and my joints are less stiff. I am so glad we met at the boat show. Thank you so much."

FREQUENTLY ASKED QUESTIONS

Is the Therapy System medically approved?

It is an approved medical device manufactured in Germany to the highest standards receiving the Berlin Certification for medical devices and the medical product conformity CE 0633.

What does the term Cyclo-ssage mean?

Cyclo-ssage is an abbreviation for Cycloidal Massage... a non-invasive multi-directional therapy action.

How does the Therapy System Work?

This portable, full body Therapy System uses a combination of 9 Cycloidal Therapy units and 6 infrared heat lenses helping to alleviate the secondary symptoms of a variety of health conditions.

How will the Cycloidal Therapy benefit me?

The non-percussive and controllable therapy action will increase blood and lymphatic circulation by stimulating deep muscle activity.



What is the likely outcome of using the Cyclo-ssage System regularly?

The combination of Cycloidal Therapy and infrared heat treatment is a very effective method of providing on-going pain and ache relief whilst also helping the body generally work more efficiently.

Are there any side effects?

There have been no findings with regard to negative secondary effects. The Cycloidal Therapy units produce a non-invasive action allowing muscles to relax and circulation to improve without aggravating any skeletal or joint problems that pre-exist.

Is it simple to operate?

A hand held digital control will enable you to time and to control the intensity of your own personal therapy treatments.

Where can the Therapy System be used?

The fully portable Therapy System can be used on any flat surface such as a bed, therapy table or even on the floor.

MORE FREQUENTLY ASKED QUESTIONS

Can I use the Therapy System daily?

Daily use is advised in order to encourage and maintain an improvement in many and varied conditions. If necessary a treatment with the Therapy System can be applied several times a day as the beneficial effects are cumulative and increase with regular treatments.

Is the Therapy System expensive to use?

No: It operates on a 12v transformed power supply using a maximum of approximately 85 watts.

Do I have to pay Vat?

The Therapy System can be zero rated for Vat if you purchase the Therapy System for your personal use and you complete a self-declaration form stating that you suffer from a qualifying chronic medical condition.



**PAIN RELIEF
MEDICATION
FREE**

What if there is a technical problem?

If you experience a technical problem with the Therapy System it is backed by a two Year Manufacturers guarantee. If you experience any difficulties or need any technical or user advice just contact our customer care department.

Do you have any customer testimonials?

You can view numerous testimonials on our website from medical professionals and customers who have found the Therapy System to be a viable alternative to prescribed medication including painkillers.

How can I find out if Cyclo-ssage can help my particular condition or needs?

We offer a personalised demonstration so that you can experience the many benefits for yourself. Please see the contact options on the back page of this brochure.

HEALTH PROFESSIONALS RECOMMEND CYCLO-SSAGE

MR PHILIP HEATON

Consultant Orthopedic Surgeon

"I have been asked to give my professional opinion with regard to the benefits gained by both my NHS and private-patients when they have been treated with the Cyclo-ssage Pro Personal Therapy System. I have written the following articles about the efficacy of using the Therapy System as part of a treatment regime for the following" ...

- Osteoarthritis
- Back Pain
- Cervical Myelopathy
- Asthma
- Osteoporosis
- Accelerated Recovery
- Thromboembolic Disease
- Anti – Ageing Therapy
- Plantar Fasciitis
- Spinal Flexibility
- Quadratus Lumborum Syndrome
- Parkinson's Disease
- Fibromyalgia
- Multiple Sclerosis
- Chronic Fatigue Syndrome
- Patient Care
- Shoulder Stiffness and Instability
- Diabetes Mellitus
- Oedema
- Stress

These articles can be viewed on the Cyclo-ssage website: www.cyclo-ssage.com

DR KATHRYN LEWIS GP

Westbury Surgery Shrewsbury

"The use of the Cyclo-ssage Pro Personal Therapy has been well received by the practices population. It has been used for a variety of medical conditions and has improved pain relief and recovery in acute conditions. It has also reduced the need for strong analgesics resulting in a reduction in prescribing costs and in patients reporting side effects."



HOW CAN I ARRANGE A CYCLO-SSAGE DEMONSTRATION?



VISIT US AT SHOWS AND EXHIBITIONS

Cyclo-ssage exhibit at over 80 shows and various other venues throughout the UK and Ireland each year. Visit our website and check our show diary to find a venue near you that we will be attending.



REQUEST A HOME DEMONSTRATION

Experience yourself how Cyclo-ssage can help your condition – go to our contact details on the back page of this brochure and we will be happy to arrange a personal demonstration for you in the privacy of your home or workplace.



WORKING WITH CYCLO-SSAGE AS AN AGENT OR A DISTRIBUTOR

If you are interested in working with Cyclo-ssage as an agent in the UK, or as an agent or distributor overseas please contact us for further information.

CYCLO-SSAGE OVERSEAS



Cyclo-ssage have distributors and agents in various countries including America, New Zealand and Norway. Please enquire to see if we have a distributor in your area and we will contact them for you.



FOR CUSTOMER CARE, FURTHER INFORMATION
OR TO MAKE AN ENQUIRY...

Email



enquiries@cyclo-ssage.com

Freephone



0800 169 2808

Website



www.cyclo-ssage.com

Social Media



Facebook & Twitter

Video



You Tube

QR Reader



Scan

DEMONSTRATIONS AVAILABLE ON REQUEST

Cyclo-ssage Ltd, Units 18-19, Sycamore Close, Duffryn Business Park, Ystrad Mynach,
Mid Glamorgan, South Wales, CF82 7RJ, U.K.

Freephone: 0800 1692808 | Office Tel: 01443 815496 | Fax: 01443 862227

Email: enquiries@cyclo-ssage.com | Vat No: 938 2636 95 | Company Reg: 6301170